




BOVIE[®] BANTAM|PRO[™]

Electrosurgical Generator

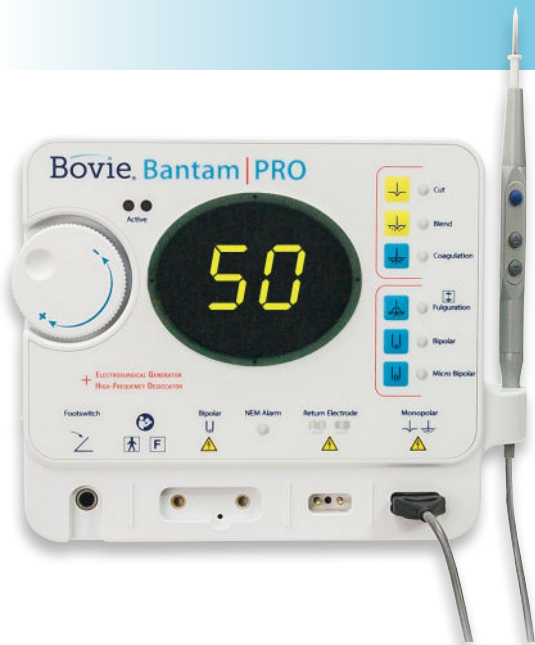

The Bovie Bantam|PRO, A952, is a High Frequency Desiccator with cutting ability. You might call it a split personality. One half of the unit is a high frequency desiccator, an upgrade to our very popular Aaron[®] 940. The other half is a Bovie Generator, giving you 50 watts of cutting power and coag capabilities – plenty of power for most cutting procedure performed in the doctor's office.

High Frequency Desiccator with 50 watts of Cutting capabilities

Two machines in one convenient package

Adjustable in 1/10th watt increments in the first 10 watts, for the most accurate output available

No grounding pad required when working in the desiccation mode

When in electrosurgical generator mode, the unit features Bovie NEM[™] pad sensing technology which monitors the return electrode for optimum patient protection.



BOVIE® BANTAM|PRO™ A952 ELECTROSURGICAL GENERATOR



STANDARD EQUIPMENT

A952 (1) electrosurgical generator,
A902 (1) handpiece
A910 (2) non-sterile disposable handpiece sheaths*
A910ST (2) sterile disposable handpiece sheaths*
ES01 (3), **ES02** (1), **ES20** (1) sterile disposable electrodes
A804 (5), **A806** (5), **A806DE** (2) disposable non-sterile dermal tips*
A805 (2), **A807** (2), **A807DE** (2) disposable sterile dermal tips*
A1252C (1) reusable grounding cord
ESRE (5) disposable split grounding pads
A837 wall mount kit*
 110 VAC Hospital-grade power cord*
 User's guide (MC-55-238-001)*

* Not Pictured



A902



ES01, ES02, ES20



A1252C



ESRE

COMPLEMENTARY PRODUCTS †



A813
Table Top Stand †



A803
Footswitch †




A827V
Forcep Cable †



A824
Reusable Bipolar Forceps †



Lletz Loop
Electrodes †

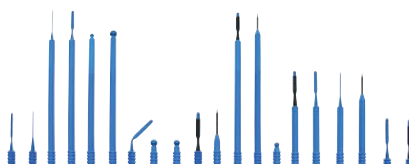

Resistick II
Coated Electrodes †



Tungsten Wire
Loop Electrodes †

† Sold Separately

PRODUCT SPECIFICATIONS

DIMENSIONS

Height: 188mm (7.4 in.)
 Width: 228mm (8.98 in.)
 Depth: 105mm (4.13 in.)
 Weight: < 2.26 kg (< 5 lbs.)

OUTPUT POWER

Cut: 50 watts
 Blend: 50 watts
 Coagulation: 50 watts
 Fulguration: 40 watts
 Bipolar: 40 watts
 Micro Bipolar: 40 watts

OUTPUT FREQUENCY

304 – 451kHz

LINE VOLTAGE

A952: 100-240 VAC

LINE FREQUENCY

50Hz

WARRANTY

Four Years