

burton

Luxo Medical Products

# OUTPATIENT<sup>®</sup> II & COOLSPOT<sup>™</sup> II

*Individually or in combination, Burton's Outpatient<sup>®</sup> II and CoolSpot<sup>™</sup> II are the right lights for minor surgery and diagnostics!*

PERFORMANCE RELIABILITY VALUE



Burton's Outpatient<sup>®</sup> II and CoolSpot<sup>™</sup> II are the perfect solutions for hospital specialists, clinics and exam rooms where versatile, high-quality illumination is required.

*The Right Light*

DIAGNOSTIC & MINOR SURGERY

# OUTPATIENT® II

*The perfect blend of minimal heat output and high color temperature, low-maintenance and top performance.*

The Outpatient® II is the documented industry leader for diagnostic and procedure lights. More physicians choose the Outpatient® II than any other diagnostic light for focusable flood-type light, maximum intensity and shadow-free operation. With minimal heat output and high color temperature, the Outpatient® II is the perfect low-maintenance, top performing light for nearly every application in the hospital, surgical suite or physician's office.



## OUTPATIENT® II FEATURES:

- > 8000 footcandles (86,000 Lux) at 24 inches (61 cm)
- > 3300° K color temperature
- > Focuses with central removable, autoclavable SteriHandle™
- > Friction knob allows user to position self-balancing arm without a tool for drift-free positioning
- > Cool operation assured with heat filters and dichroic coatings
- > Shadow-free operation provided by 3 "optically improved" 50-watt halogen bulbs, each with 5600 hr. average bulb life
- > Bulbs are wired in parallel to assure continuous operation
- > Design allows removal or replacement of SteriHandle™ with one hand (for assuring sterile procedures)
- > Floorstand units include footswitch, locking casters and cord wrap
- > Mounting systems meet California Seismic Codes
- > 115V model standard; 230V models also available
- > 5-year warranty
- > UL 60601-1/60601-2-41 certified
- > Made in the USA

## MODELS (complete models only)

115V, 60Hz	230V, 50Hz	
#OP216FL	#O202FL	Outpatient® II Floorstand
#OP216SC	#O202SC	Outpatient® II Single Ceiling
#OP216DC	#O202DC	Outpatient® II Double Ceiling
#OP216W	#O202W	Outpatient® II Wall
#OP216ST		Single Outpatient® II on Fastrac™
#OP216DT		Double Outpatient® II on Fastrac™

## ACCESSORIES

#0006130PK	Set of 3 Hi-intensity 50-watt Replacement Bulbs
#1008428	Outpatient® II SteriHandle™
#0008100PK	Set of 25 Burton Sterile Disposable Light Handle Covers



# COOLSPOT™ II

*The advanced health care lighting solution that fits  
both your space and your budget.*



The aesthetically appealing CoolSpot™ II offers incredibly pure, white light along with a guaranteed no-drift arm, high intensity setting and variable spot size control. The CoolSpot II offers intensity of up to 13,000 footcandles at 24 inches (61 cm) and 5000 footcandles at 1 meter. The CoolSpot™ II Minor Surgery Light is the spotlight of choice.

## MODELS (complete models only)

115V, 60Hz	230V, 50Hz	
#CS316FL	#C302FL	CoolSpot™ II Floorstand
#CS316SC	#C302SC	CoolSpot™ II Single Ceiling
#CS316DC	#C302DC	CoolSpot™ II Double Ceiling
#CS316W	#C302W	CoolSpot™ II Wall
#CS316ST		Single CoolSpot™ II on Fastrac™
#CS316DT		Double CoolSpot™ II on Fastrac™

## ACCESSORIES:


#0007006PK	Set of 4 EKE Replacement Bulbs
#1007696	Universal SteriHandle™ Replacement Model
#0008100PK	Set of 25 Burton Sterile Disposable Light Handle Covers

## COOLSPOT™ II FEATURES:


- > Up to 13,000 footcandles (140,000 Lux) at 24 inches (61 cm) and 5000 footcandles (54,000 Lux) at 1 meter
- > Up to 4300° K color temperature
- > Variable intensity settings and high intensity switch
- > Variable spot size of 2"-10" (5-25cm)
- > New, improved and easy to clean housing design with convenient SteriHandle™ for positioning
- > Removable, autoclavable center SteriHandle™ design allows removal or replacement of SteriHandle™ with one hand (for assuring sterile procedures)
- > 150-watt halogen bulb provides 400 hr. average bulb life on standard setting
- > Friction knob allows user to position self-balancing arm without a tool for drift-free positioning
- > Floorstand units include footswitch, locking casters and cord wrap
- > Mounting systems meet California Seismic Codes
- > 115V model standard; 230V models also available
- > 5-year warranty
- > UL 60601-1/60601-2-41 certified
- > Made in the USA



# Compare technical data and you will choose Burton's OUTPATIENT® II & COOLSPOT™ II

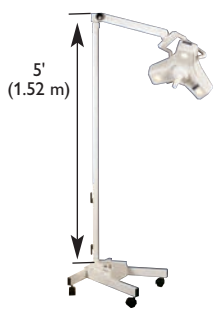


#OP216FL Floorstand  
Now, includes footswitch,  
locking casters & cord wrap  
(See Notes 1, 4, 5)

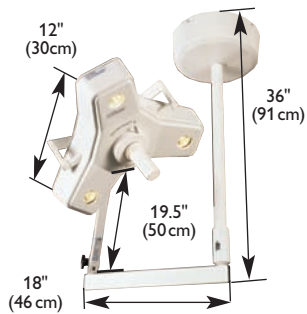


#OP216SC Single Ceiling  
(See Notes 1, 2, 3, 4, 5)




#OP216W Wall  
(See Notes 1, 4, 5)

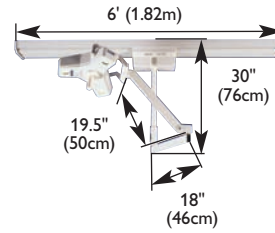

#OP216ST Single Fastrac™  
(See Notes 1, 2, 3, 4, 5)

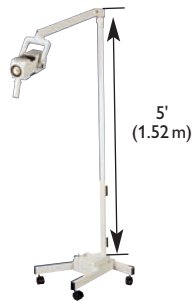


#CS316FL Floorstand (See Notes 1, 4, 5)  
Now includes footswitch,  
locking casters & cord wrap

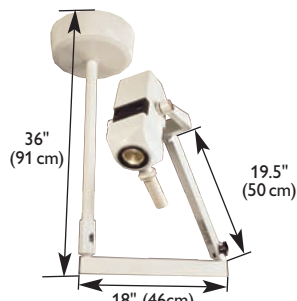


#CS316SC Single Ceiling  
(See Notes 1, 2, 3, 4, 5)




#CS316W Wall Mount  
(See Notes 1, 4, 5)

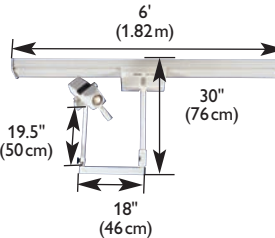

#CS316ST Single Fastrac™  
(See Notes 1, 2, 3, 4, 5)

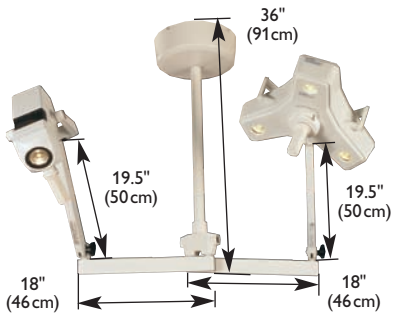

#OCS16CC CoolSpot™ II / Outpatient® II  
with Double Ceiling Mount (See Notes 1, 2, 3, 4, 5)

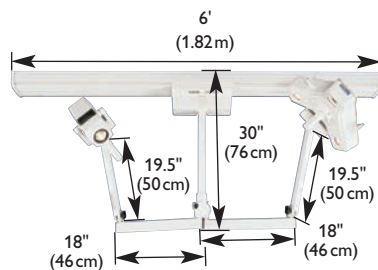

#OCS16CT CoolSpot III/Outpatient II Combo  
with Fastrac™ (See Notes 1, 2, 3, 4, 5)

Note 1: Standard dimensions in inches (centimeters) shown. All ceiling mounts shown on standard 9' (2.74 m) ceilings.

Note 2: Options available for custom lengths for both Fastrac™ and drop tubes. Ask for Burton's customer service for custom Fastrac lengths. For single/double mounts, either drop tube of 30" (76 cm) or 42" (1.07 m) can be specified at time of order only. If longer drop tube is required, order #1008060 and specify length (up to 58" (1.47 m)).

Note 3: Standard 30" (76 cm) drop tube is generally acceptable for 9' (2.74 m) ceilings. If shorter drop tube is required, Burton will manufacture (at no cost) requested length at time of order only. Special lengths not specified at time of order are subject to additional charges.

Note 4: The optical data are nominal values based on measurements done according to the UL/IEC performance standards at 120V at 24"

Note 5: Contact Burton for part #'s on 115V, 50Hz/230V, 60Hz models.



**burton**

Luxo Medical Products

BURTON MEDICAL  
21100 Lassen Street  
Chatsworth, CA 91311  
Phone: 800-444-9909, 818-701-8700  
Fax: 800-765-1770, 818-701-8725  
www.burtonmedical.com

The Right Light