

# care-e-vac<sup>®</sup>3

PORTABLE ASPIRATOR



High Flow and High Suction

## DESCRIPTION AND APPLICATION

Care-e-vac 3 is high-vacuum and high-flow AC/DC powered portable suction aspirator. The aspirator's high suction capacity and its easy handling assist medical personnel carrying out effective and proven patient procedures. It is an excellent choice for hospital code carts, patient transport and wherever strong battery powered suction is needed.

Care-e-vac 3 contains a rechargeable battery monitoring system that is designed to hold the battery at optimal charge, without harmful overcharging. The internal battery monitoring system will:

- Bring the battery to full charge
- Maintain "float" charge when the battery reaches full charge, eliminating damage to the battery caused by overcharging
- Monitor and display battery charge level during operation

## KEY FEATURES


- High suction capacity up to 50 L/min flow on AC power and 45 L/min on battery power at full vacuum level
- Rechargeable battery with long run time of more than 60 minutes
- Safe secretion collection canister with overflow protection and disposable hydrophobic filter
- Adjustable canister bracket to accommodate different size canisters
- Robust housing made of impact-resistance polycarbonate and ABS materials

## IDEAL FOR

- Pharyngeal Aspiration
- Gastric Suctioning
- Tracheal Suctioning
- Code Carts
- Patient Transport
- Nursing Homes
- Homecare
- Physician's Office
- Surgery Centers




## TECHNICAL SPECIFICATIONS

Part Number	US: 758000	International: 758010
Electrical Rating	120 VAC 0.75A, 60Hz 	100-240 VAC 0.75A, 50/60Hz
Battery	12V 5AH Rechargeable Sealed Lead Acid Battery operating time: 60 minutes (Uninterrupted with fully charged/new battery) Charging time: 6 hours to 85%, 8 hours to 100%	
Pump	Oil-less diaphragm	
Suction Capacity	Vacuum: 550 mmHg at full vacuum on AC power Flow: 50 L/min on AC power at full vacuum 45 L/min on DC power at full vacuum with fully charged battery	
Vacuum Settings	Variable with regulator on panel	
Vacuum Gauge	Range 0-760 mmHg with 3% Accuracy of full scale	
Collection Canister	800cc disposable plastic with mechanical shutoff	
Canister Bracket	Adjustable – will accept up to 1200cc standard canister	
Suction Tube	Disposable clear plastic tubing with ¼" ID and 6' long with rubber barb connector on both ends	
Suction Tip	Standard bulb tip with barb connector	
Filter	Disposable Hydrophobic filter with 14" tubing	
Manufacturer's Warranty	3 years	
Dimensions (Size & Weight)	8.5" x 10.50" x 8.5" (Without canister bracket) (10.5 lbs.)	
Shipping (Size & Weight)	11 ½" x 18 ½" x 11 ½" (15 lbs.)	

## ADDITIONAL FEATURES AND BENEFITS

- Constant suction providing 550mmHg and 50 L/min free air flow
- Easy to clean and process for quick turnaround
- Compact and lightweight
- Crash cart mount (optional)
- 3-year warranty
- Permanently attached power cord

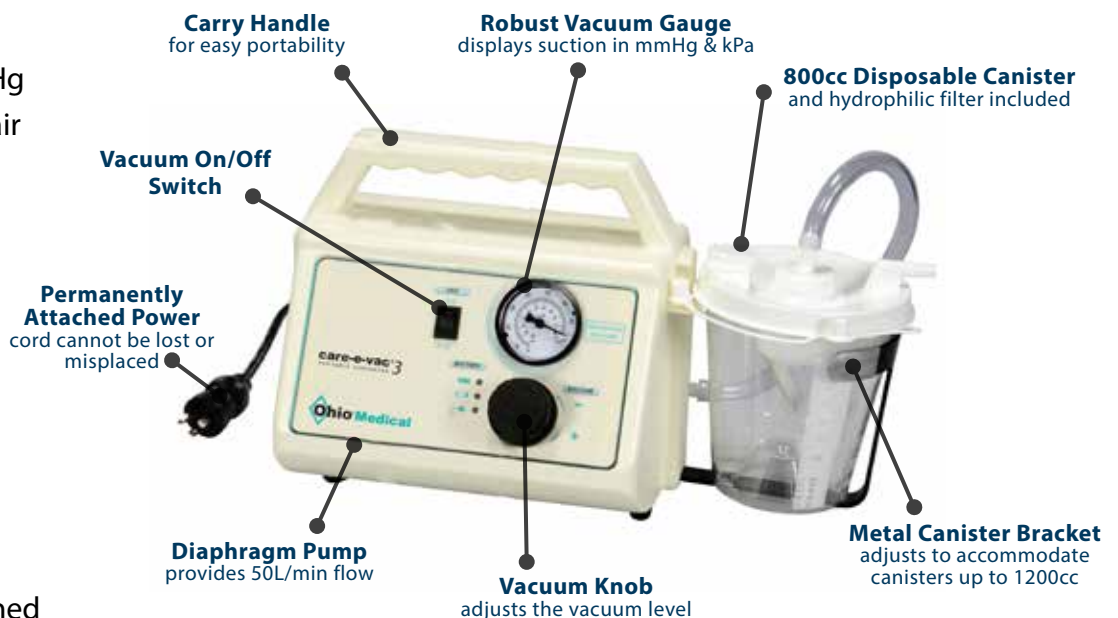

# AVAILABLE ACCESSORIES

758007	Crash Cart Mount Adapter Kit
AI2557ASY	Crash Cart Mount (requires adapter kit)
AI2580000	Floor Stand (requires adapter kit)
749144	Replacement Sealed Lead Acid Battery
AI4080000	Collection Canister, 800cc
AI0480010	Canister Only Repl. Kit, 800cc, Lid Included (Pkg of 10)
AI4080025	Canister Only Repl. Kit, 800cc, Lid Included (Pkg of 25)
AI4120000	Collection Canister, 1200cc
AI4120010	Canister Only Repl. Kit, 1200cc, Lid Included (Pkg of 10)
AI4120025	Canister Only Repl. Kit, 1200cc, Lid Included (Pkg of 25)
751102	Kit, Hydrophobic Filter/14" Tubing
749002	Tip, Yankauer Suction
749003	Suction Tubing, 72" x 1/4" ID
750810	Kit, 800cc, Aspirating Change-out (Case of 10)
750825	Kit, 800cc, Aspirating Change-out (Case of 25)

\*\*\*If you need a cord other than the standard U.S. power cord, please add the desired cord part number (Example: 758010/751042 with Italian plug)

751041 UK/Ireland	
751042 Italian	
751043 Japanese	
751044 Australian	
751045 Israel	
751046 CEE 7/7 plug	
751024 U.S. Power Cord	



© 2017 Ohio Medical LLC. All rights reserved. This document contains information that is proprietary and confidential to Ohio Medical LLC. Use of this information is under license from Ohio Medical LLC. Any use other than that authorized by Ohio Medical LLC is prohibited. Ohio Medical and the Ohio Medical logo are registered trademarks of Ohio Medical LLC.

04/17 (Rev. 6)  
255391