

 symmetry surgical® | **Bovie**®
HIGH FREQUENCY DESICCATORS



BOVIE® DERM 942 HIGH FREQUENCY DESICCATOR




The DERM 942 and its patent-protected handpiece are designed and engineered to be the most reliable and durable desiccator available today. The unit's life is prolonged by incorporating the latest in chip technology and reducing the number of internal parts, thereby generating less internal heat.

The DERM 942 comes with 40 watts of coagulation capabilities and is adjustable in 1/10th watt increments in the first 10 watts (in monopolar and bipolar). It is a ground reference unit which means no grounding pad is required. The DERM 942 is the most cost-effective and user-friendly unit on the market, meaning less to do and less to buy. Set-up and usage are very simple and straightforward when compared to the competitor. For example, adjusting from low to high power and adding foot control to the handpiece are accomplished in the following ways.


Bovie DERM 942	Competitor
Adjusting from Lo to Hi Power	
1-Step Adjust rotary dial or press handpiece buttons	4-Steps Unplug handpiece from Lo power receptacle Move power selector switch from Lo to Hi Plug handpiece into Hi power receptacle Adjust rotary dial or press handpiece buttons
Adding Foot Control to Handpiece	
1-Step Plug footswitch into front of unit	6-Steps Unplug handpiece from bottom of unit Unplug handpiece from front of unit Detach handpiece from holder Attach foot control handpiece to holder Plug foot control handpiece into front of unit Plug footswitch into bottom of unit

BOVIE® DERM 942 HIGH FREQUENCY DESICCATOR

PRODUCTS INCLUDED

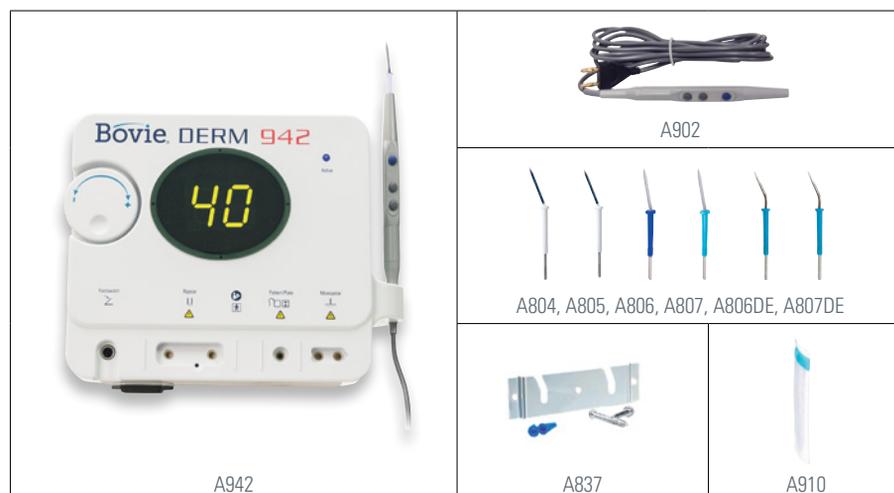

▼ COMPLEMENTARY PRODUCTS

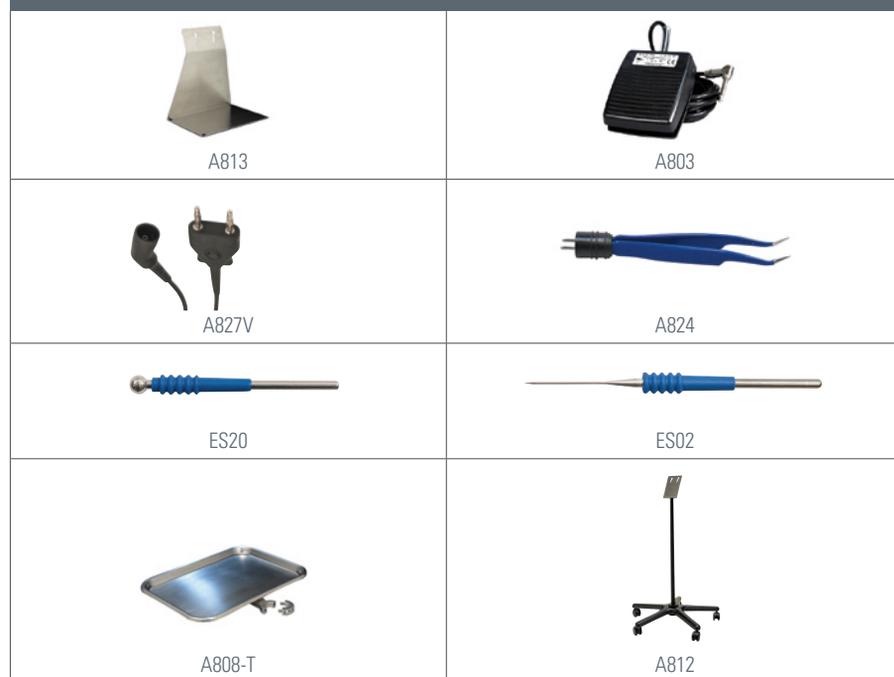

Style # A942

▼ Products Included

- A942** (1) high frequency desiccator
- A902** (1) handpiece
- A804** (5) non-sterile sharp dermal tips
- A805** (2) sterile sharp dermal tips
- A806** (5) non-sterile blunt dermal tips
- A807** (2) sterile blunt dermal tips
- A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips
- A807DE** (2) sterile DERM-Elite Premium blunt dermal tips
- A910** (2) non-sterile disposable handpiece sheaths
- A837** (1) wall mount kit, 110 VAC Hospital-grade power cord*
User's guide (**MC-55-239-001**)*

Dimensions

Height: 188mm (7.40")
Width: 228mm (8.98")
Depth: 105mm (4.13")
Weight: < 2.26kg (< 5lbs.)

Output Power

To 40 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

Line Voltage

100-240 VAC

Output Frequency: 328-492 kHz

Line Frequency: 50-60 Hz

Warranty: FOUR YEARS

*Not Pictured

†Sold separately

BOVIE® BANTAM | PRO ELECTROSURGICAL GENERATOR



The Bovie® Bantam PRO A952 is a high frequency desiccator that also has cutting ability. The unit provides 50 watts of cutting power and coagulation modes – plenty of power for most procedures performed in the doctor's office. Power is adjustable in 1/10th watt increments from 0 to 10 watts, for the most accurate output available. No grounding pad is required when working in desiccation mode.

The unit also features NEM™ (Neutral Electrode Monitoring) pad sensing technology which monitors the return electrode for optimum patient protection. The Bovie® Bantam PRO A952 is backed by Bovie's four (4) year manufacturer's warranty.

BOVIE® BANTAM | PRO ELECTROSURGICAL GENERATOR

PRODUCTS INCLUDED

 <p>A952</p>	
 <p>A902</p>	 <p>A804, A805, A806, A807, A806DE, A807DE</p>
 <p>ES01, ES02, ES20</p>	 <p>A837</p>
 <p>A910</p>	 <p>A910ST</p>

Style # A952

▼ Products Included

- A952** (1) electrosurgical generator
 - A902** (1) handpiece
 - A804** (5) non-sterile sharp dermal tips
 - A805** (2) sterile sharp dermal tips
 - A806** (5) non-sterile blunt dermal tips
 - A807** (2) sterile blunt dermal tips
 - A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips
 - A807DE** (2) sterile DERM-Elite Premium blunt dermal tips
 - ES01** (3) blade sterile disposable electrodes
 - ES02** (1) needle sterile disposable electrodes
 - ES20** (1) ball sterile disposable electrodes
 - A910** (2) non-sterile disposable handpiece sheaths
 - A910ST** (2) sterile disposable handpiece sheaths
 - A1252C** (2) reusable grounding cord*
 - ESRE** (5) disposable split grounding pads*
 - A837** (1) wall mount kit
- 110 VAC Hospital-grade power cord*
User's guide (MC-55-238-001)*

Dimensions

Height: 188 mm (7.40 in.)
Width: 228 mm (8.98 in.)
Depth: 105 mm (4.13 in.)
Weight: < 2.26 kg (< 5 lbs.)

Output Power

Cut: 50 watts
Blend: 50 watts
Coagulation: 50 watts
Fulguration: 40 watts
Bipolar: 40 watts
Micro Bipolar: 40 watts

Output Frequency: 304-451 kHz

Line Frequency: 50-60 Hz

Line Voltage: 100-240 VAC

Warranty: **FOUR YEARS**

BOVIE® BANTAM | PRO A952-G ELECTROSURGICAL GENERATOR



The Bovie® Bantam PRO | A952-G is packaged with the essential accessories needed to perform an OB/GYN LEEP procedure. The full system is situated on a mobile stand for effortless portability. It is an extremely affordable electrosurgical generator package with smoke evacuation and includes a variety of electrodes and LLETZ loop electrodes. The unit is backed by Bovie's four (4) year manufacturer's warranty.

BOVIE® BANTAM | PRO A952-G ELECTROSURGICAL GENERATOR

PRODUCTS INCLUDED

 <p>A952</p>			
 <p>A902</p>	 <p>A1252C</p>	 <p>SE02</p>	
 <p>A804, A805, A806, A807, A806DE, A807DE</p>		 <p>A910ST</p>	 <p>A812C</p>
 <p>ES01, ES02, ES20</p>	 <p>ESLK</p>	 <p>ES07</p>	
 <p>786T</p>	 <p>A837</p>	 <p>ESRE</p>	 <p>A910</p>
 <p>SERF</p>			

Style # A952-G

▼ Products Included

▼ Electrosurgery Equipment

- A952** (1) electrosurgical generator
 - A902** (1) handpiece
 - A804** (5) non-sterile sharp dermal tips
 - A805** (2) sterile sharp dermal tips
 - A806** (5) non-sterile blunt dermal tips
 - A807** (2) sterile blunt dermal tips
 - A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips
 - A807DE** (2) sterile DERM-Elite™ Premium blunt dermal tips
 - ES01** (3) blade sterile disposable electrodes
 - ES02** (1) needle sterile disposable electrodes
 - ES20** (1) ball sterile disposable electrodes
 - A910ST** (2) sterile disposable handpiece sheaths
 - A910** (2) non-sterile disposable handpiece sheaths
 - A1252C** (2) reusable grounding cord
 - ESRE** (5) disposable split grounding pads
 - A837** (1) wall mount kit
 - ES07** (5) extended ball electrodes
- 110 VAC Hospital-grade power cord*
User's guide (**MC-55-238-001**)*

ESLK electrosurgical loop kit

- ES09** (1) 10mm x 10mm loop
- ES11** (1) 15mm x 10mm loop
- ES12** (1) 20mm x 8mm loop
- ES13** (1) 20mm x 15mm loop
- ES16** (1) 10mm x 8mm square


A812-C Mobile stand

- A812** (1) mobile stand for A952
- A812-E** (1) height extender for A812
- A808-T** (1) top tray and clamp for A812
- A812-BT** (1) bottom tray for A812

Smoke evacuation system

- SE02** (1) Smoke Shark® II with pneumatic footswitch
 - SF35** (1) 35-hour filter
 - 786T** (1) 7/8" (22.2mm) tube x 6' (182.9cm) long
 - SERF** (1) reducer fitting (attaches to vaginal speculum)
- 110 VAC Hospital-grade power cord*
User's guide (**MC-23088**)*

BOVIE® DERM 102 HIGH FREQUENCY DESICCATOR

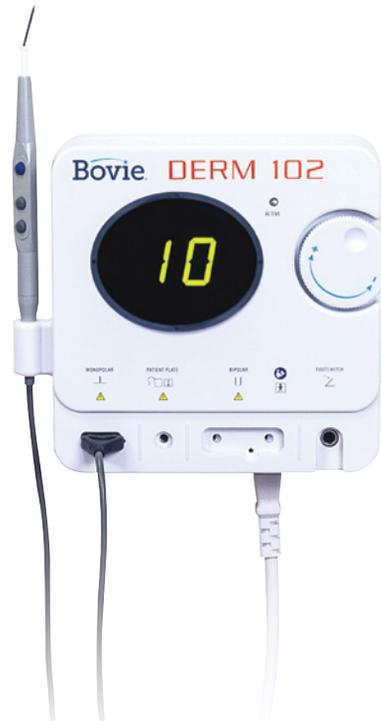

The Bovie® DERM 102 High Frequency Desiccator offers ten watts of monopolar power, adjustable in 1/10th watt increments, to desiccate and fulgurate tissue. The power can be adjusted using the 3-Button Power Control Handpiece A902 (Patent #8,998,899), or using the rotary dial on the unit. The unit features a large, bright, LED display for excellent visibility across the procedure room. If foot control activation is preferred while using the DERM 102, just plug in the A803 Bovie Footswitch (sold separately) and activate. No separate handpiece is required.

For the practitioner who wants to also use bi-polar electrosurgery along with high frequency desiccation, the DERM 102 High Frequency Desiccator provides 10 watts of bipolar power, adjustable in 1/10th watt increments. To utilize this feature, the A827V* bipolar cable, a bipolar forcep* and footswitch activation are required.

The DERM 102 is packaged ready to use. The unit comes with a total of 18 disposable dermal tips. By the time the practitioners utilizes all 18 electrodes, they will have paid for the unit (twice), and generated income for their practice.

BOVIE® DERM 102 HIGH FREQUENCY DESICCATOR

PRODUCTS INCLUDED



Style # DERM 102

▼ Products Included

- DERM 102** (1) high frequency desiccator
A902 (1) handpiece
A910 (2) non-sterile, disposable handpiece sheaths
A804 (5) non-sterile single-use sharp dermal tips
A805 (2) sterile single-use sharp dermal tips
A806 (5) non-sterile single-use blunt dermal tips
A807 (2) sterile single-use blunt dermal tips
A806DE (2) non-sterile DERM-Elite™ Premium blunt dermal tips
A807DE (2) sterile DERM-Elite Premium blunt dermal tips
A837 (1) wall mount kit
 110 VAC Hospital-grade power cord*
 User's guide CD (MC-55-229-001)*

Dimensions

Height: 208mm (8.2")
 Width: 184mm (7.2")
 Depth: 69mm (2.7")
 Weight: < 1.4kg (< 3 lbs.)

Output Power

Up to 10 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

Line Voltage

100-240 VAC

Line Frequency: 50 – 60 Hz

Main Current: 320 mA Max

Warranty: FOUR YEARS

* Not Pictured



3034 Owen Drive | Antioch, TN 37013 | 800.251.3000 (US) | +41 (0) 52 557 59 00
customerservice@symmetryurgical.com | [symmetryurgical.com](https://www.symmetryurgical.com)
©2022 Symmetry Surgical Inc. All rights reserved.
Bovie High Frequency Dessicator Brochure 03.23.2022 v1